

1 cm Thick Silicon(Ge)Detectors at Liquid Helium Temperature

From Surface to Volume (Neutrino Physics)

C.Braggio¹, G.Bressi², M.Boscardin³


G.Carugno¹, G.Galeazzi⁴, N.Zorzi³

1) INFN Sez. Padova

2) INFN Sez. Pavia , 3) IRST-ITC Trento

4) INFN Lab.Naz.Legnaro


1.Large Drifting Distances: Requirements

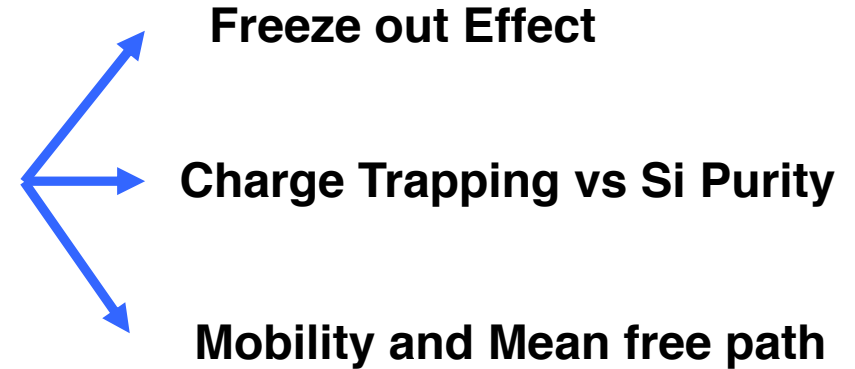

2.Silicon Processing Stages


3.Experimental Set-Up

4.Results:

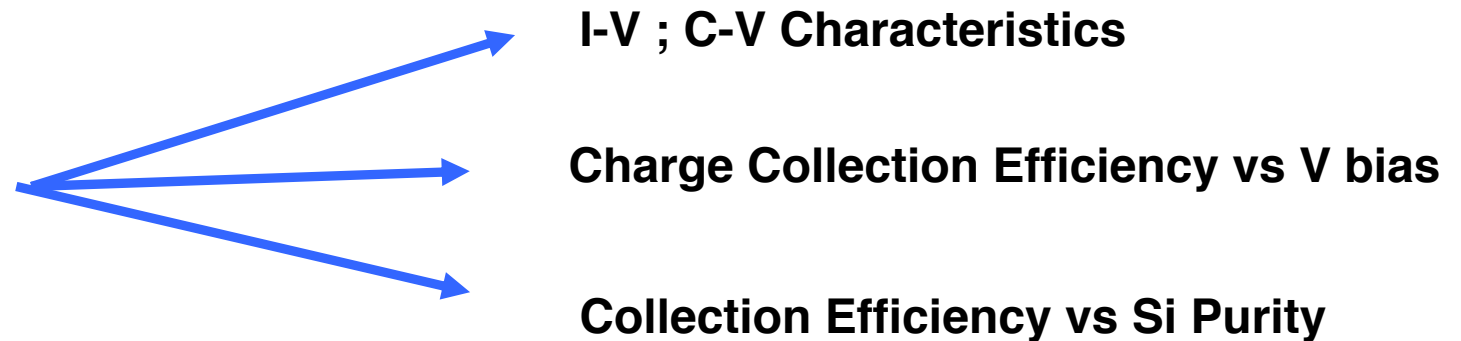

5.Conclusions

Freeze-out Effect at Low Temperature

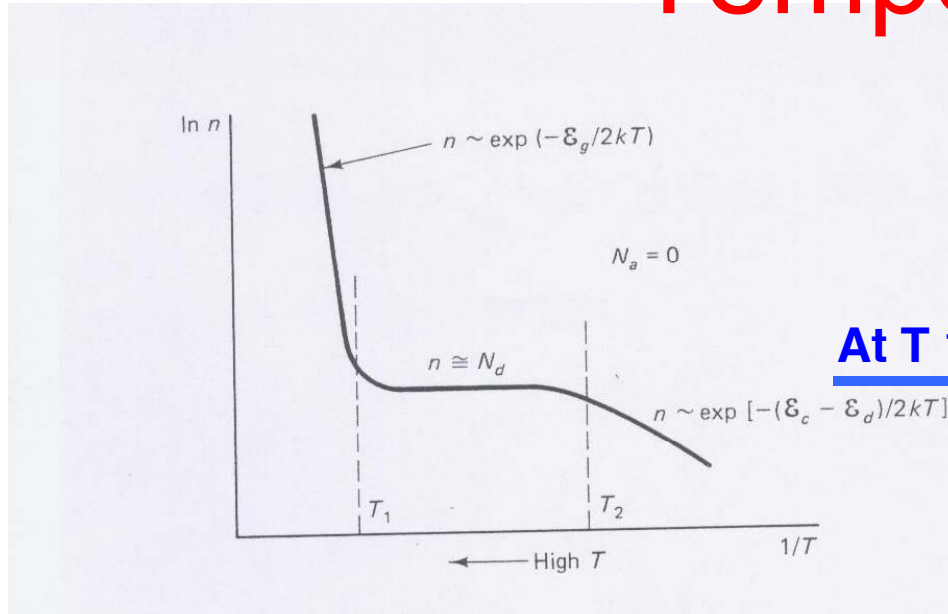

$$E_c - E_d = 10 \text{ meV} , E_a - E_v = 30 \text{ meV}$$

$$E_{KT} @ 4 \text{ Kelvin} = 0.3 \text{ meV}$$

At T 10 kelvin Silicon is Intrinsically depleted


Cap. Detector = Geometrical Cap.

BENEFITS:


No leakage current

Low Bias Voltage

But not enough

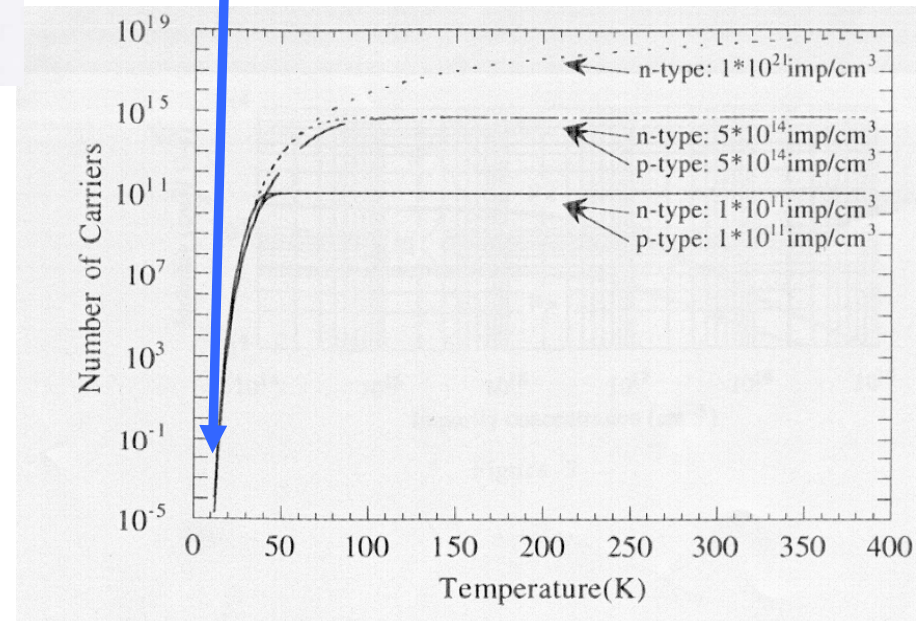

Charge Collection + Trapping at Low Temperature

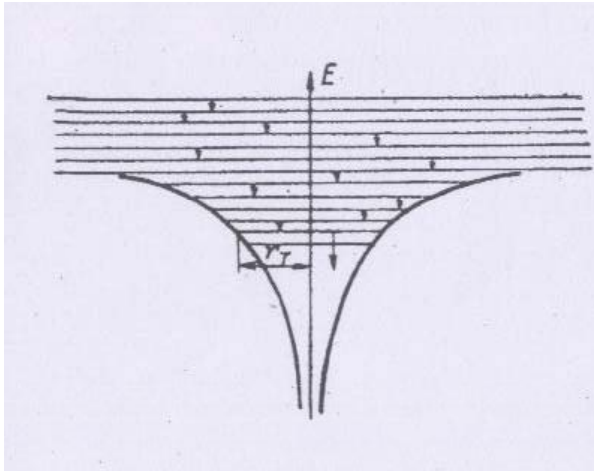


Cascade Capture Model Developed by Pitaevskii
Review Paper from Abakumov Sov.Phy.Sem.1978

2 Type of competing processes:

- a) Scattering Cross Section from Coulomb center
- b) S.C.S from “forming dipoles”

For High Purity Material (less 10^{14} imp/cm³) process b is dominating

b) Type Cross Section is proportional to $(N_d, a)^{3/4}$

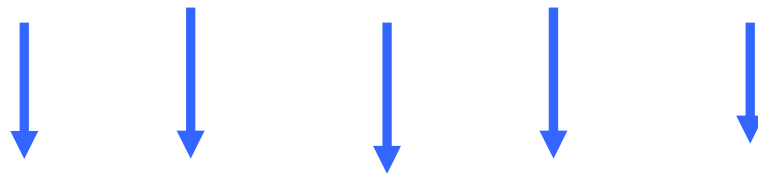

To Decrease Dipoles SCS We Need High Purity Silicon

Silicon Detector study at low Temperature for Dark Matter



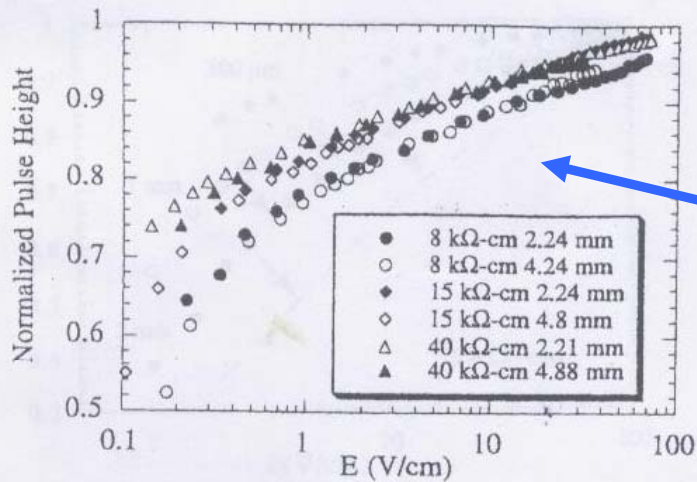


FIG. 3. Normalized pulse height vs electric field for high-purity samples.

Schotky type junctions @ 30 mKelvin


Cabrera et J.Appl.Phys.79-1996-8179

N-type Pin diode 1 mm thick resist.=10KOhmcm

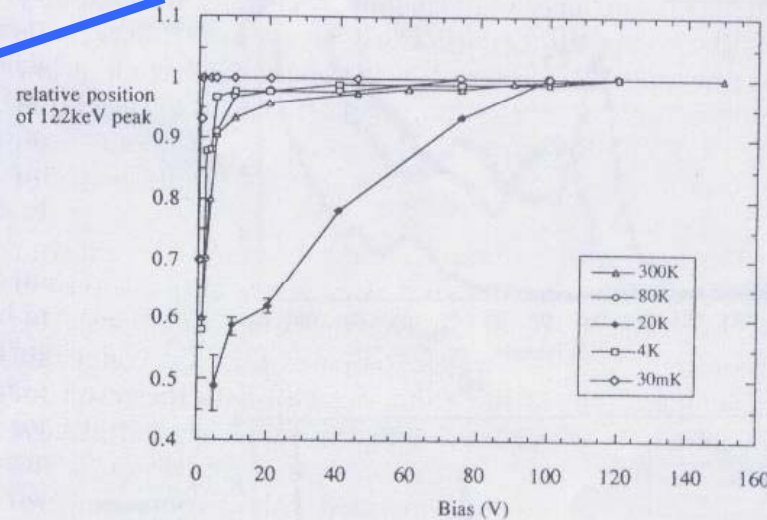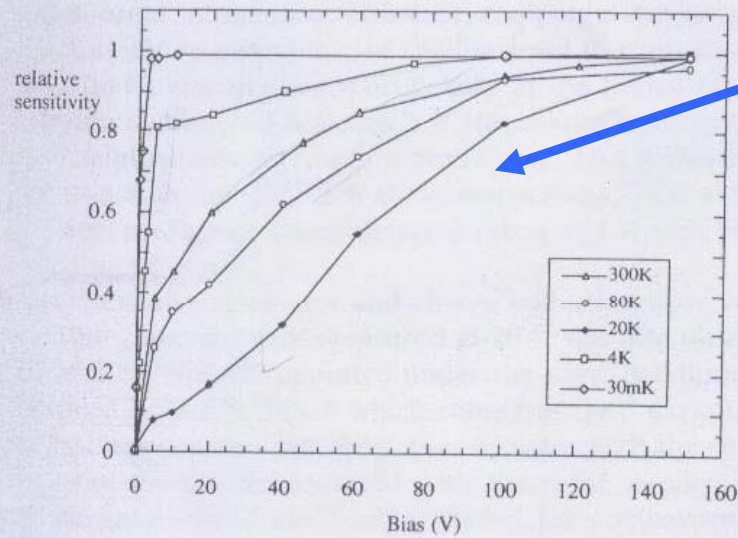

Spooner et al IEEE Trans.Nucl.Sci.40-1993-105

Mobility and Mean Free Path

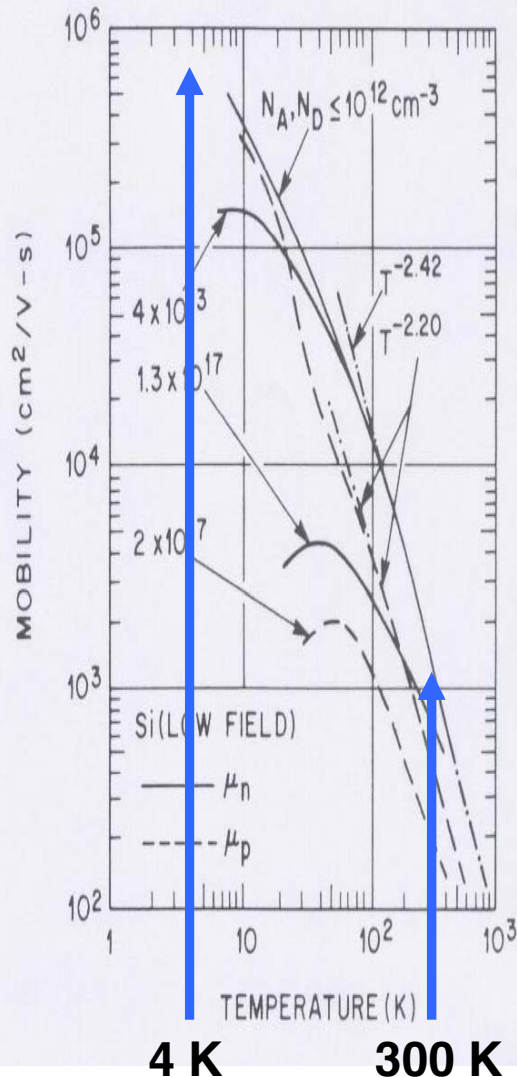

At LHe Temp. mobility increases almost 3 O.M. respect to R.T.

Electron mean free path increases at Low Temp.

$q E (\text{V/cm})$ $I = E \text{ ion (eV)}$

$I = \text{electron mean free path}$


$I = V_{\text{sat.}} * \tau$

$\tau = \text{mobility} * \text{mass} / e$

Mobility increases more than a factor 100

E avalanche from 10^5 to 10^3 V/cm

Silicon Processing Stages


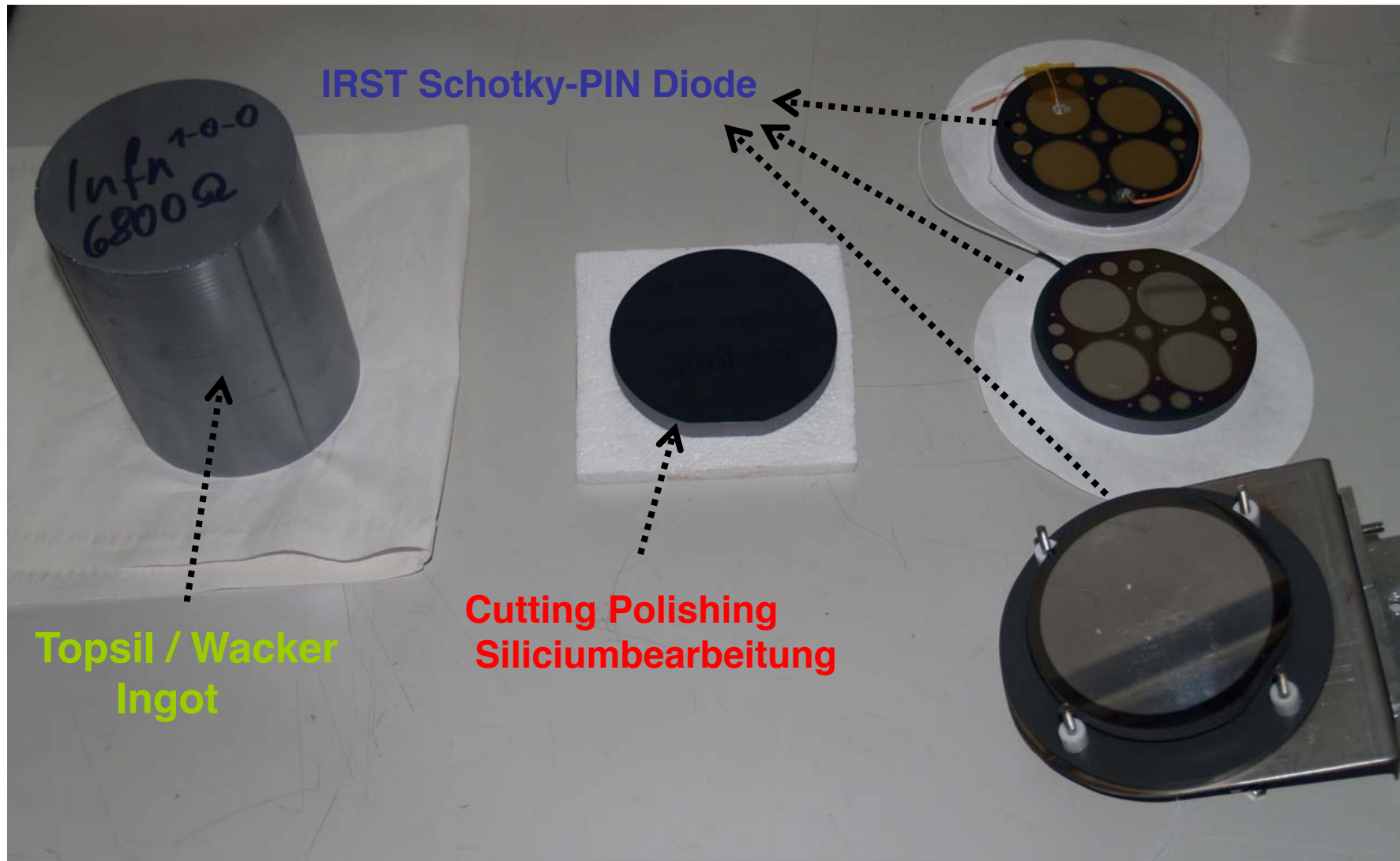


Experimental Set-Up

20 Lt LHe Dewar with LN2 Shielding

Liq. Helium Evaporation Rate: 1,3 Lq.Liter/hour




Silicon detector chamber
Filled up with 2 mbar He Gas
To attain good thermalization

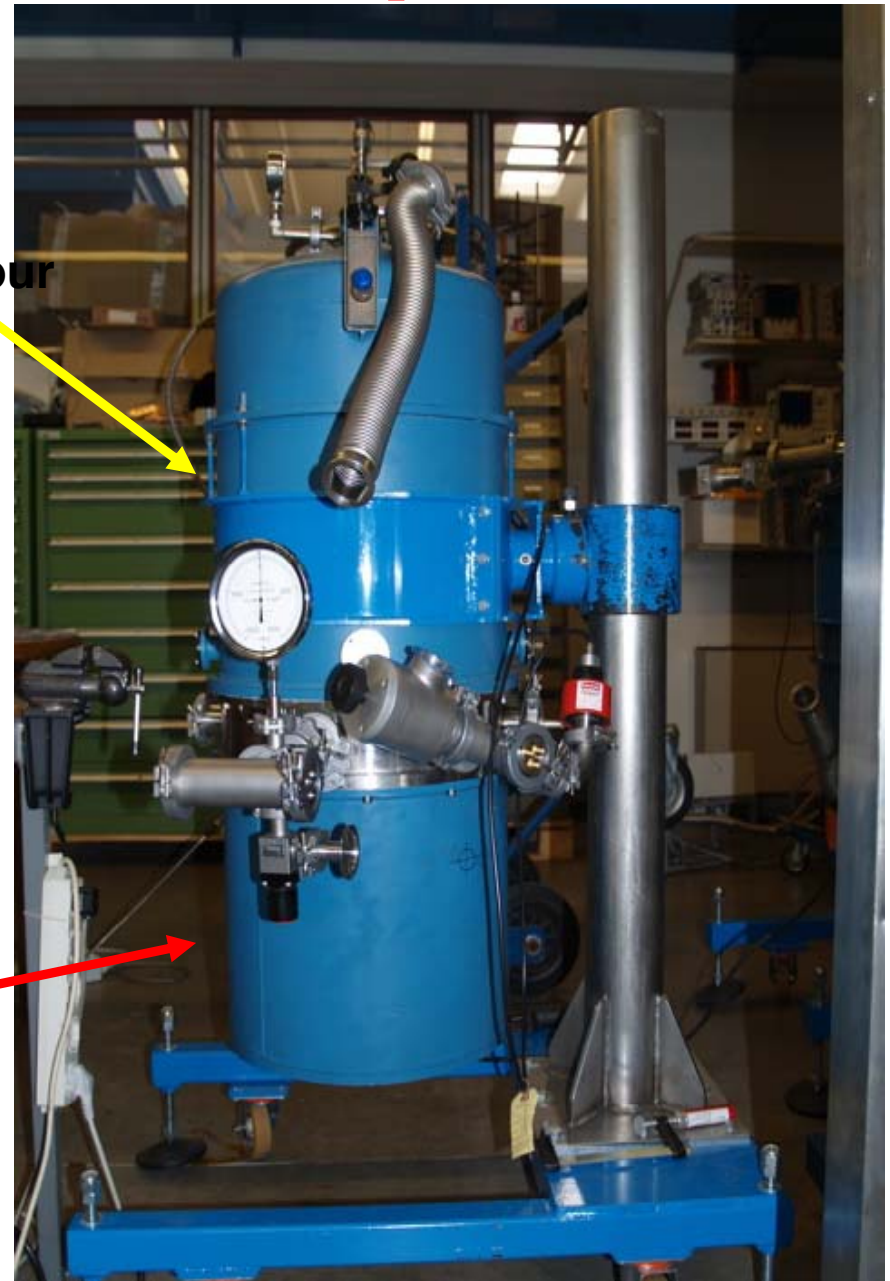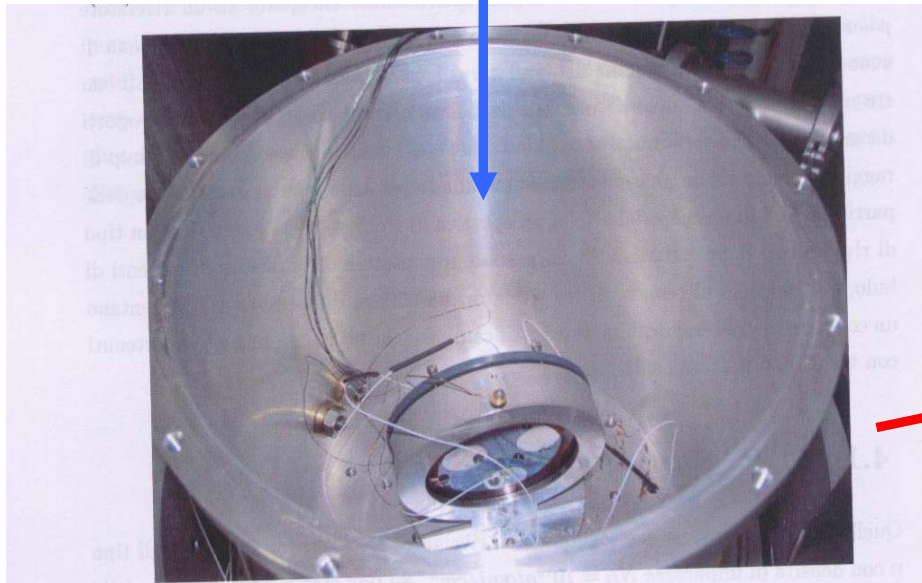


P-Type Silicon PiN diode

Preliminar meas on test diodes
and single PAD diodes

Topsil material 40 KOhm*cm RT

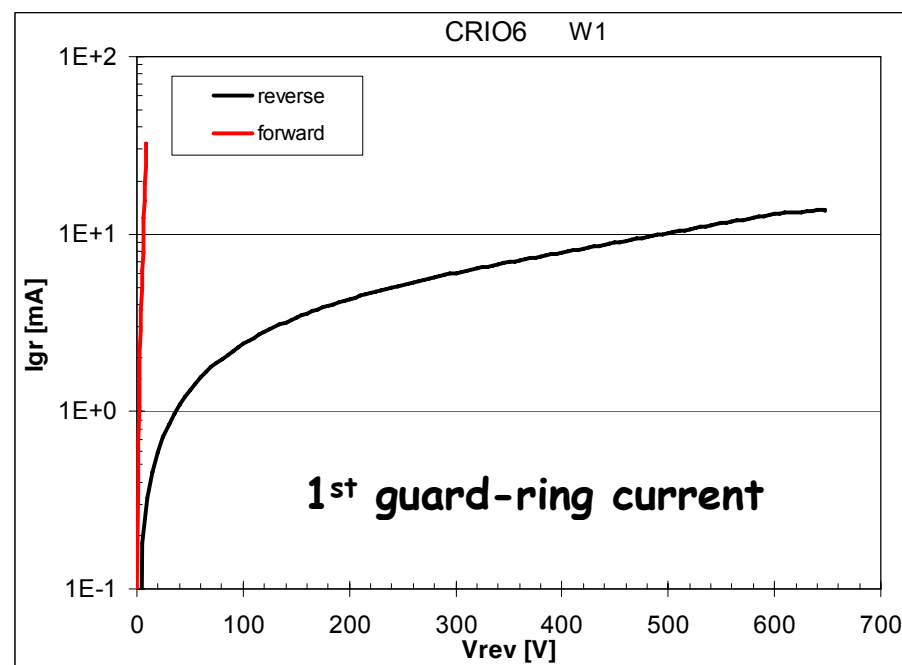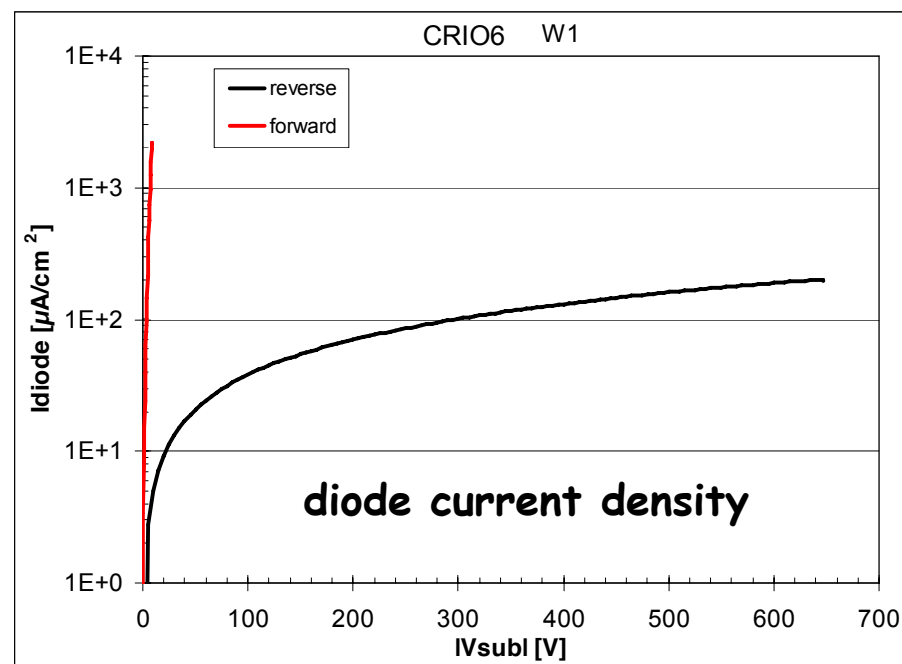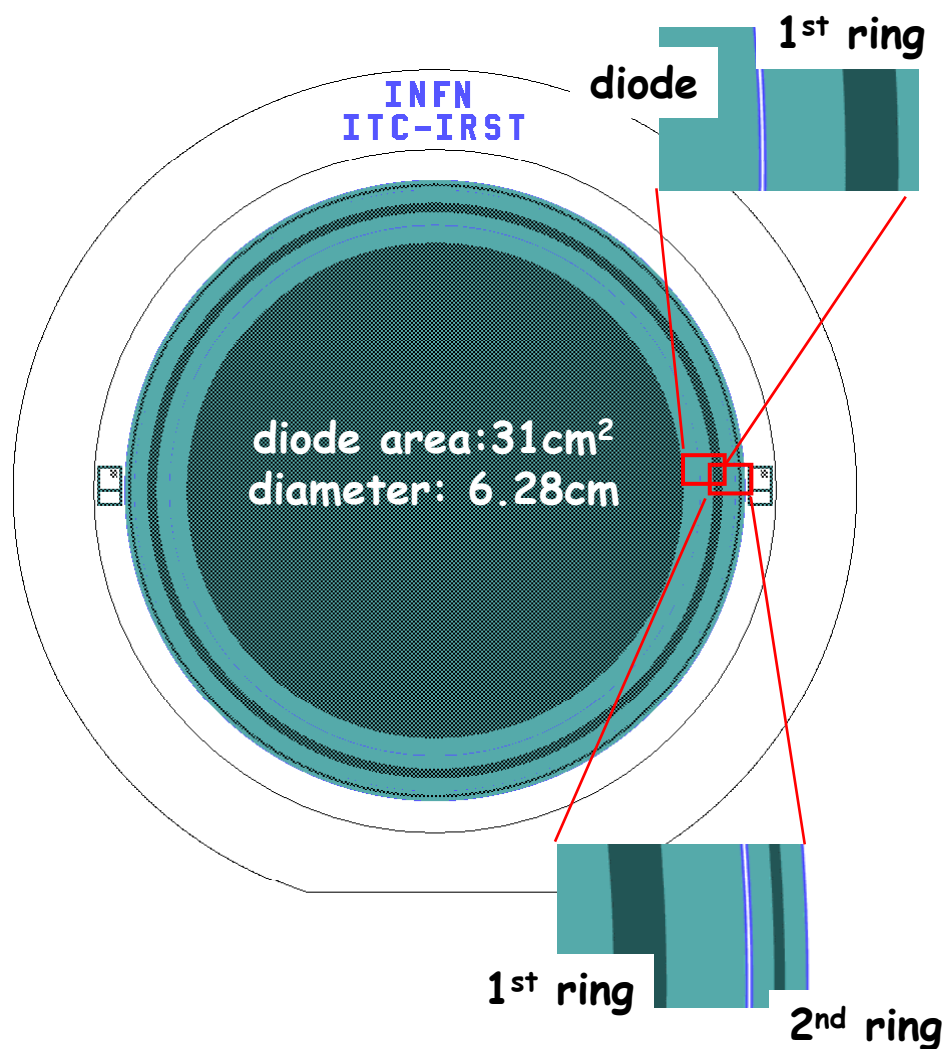

Capacitance + Leakage Current

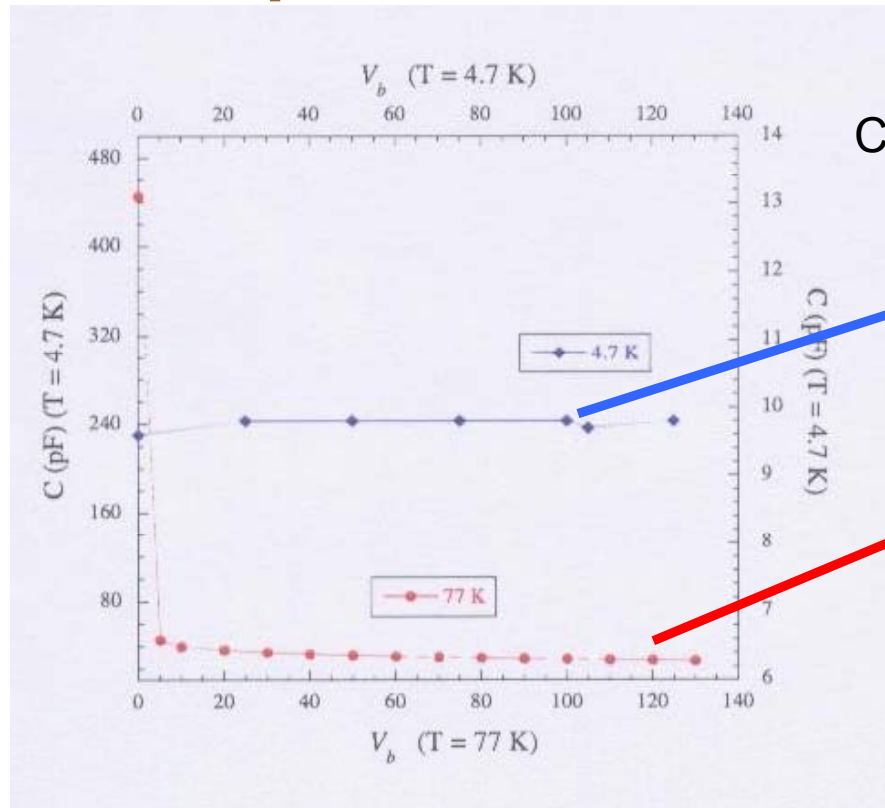


C and I vs V_{bias} measurements for 3 cm diam.


Capacity is constant vs V_{bias} at 4 Kelvin
 $C = 9.5 \text{ pF} = \text{Geom. Cap.}$

At 77 K C doesn't saturate
 @120V $C = 30 \text{ pF}$

Detector Leakage Current 1 pA/cm² at 100V

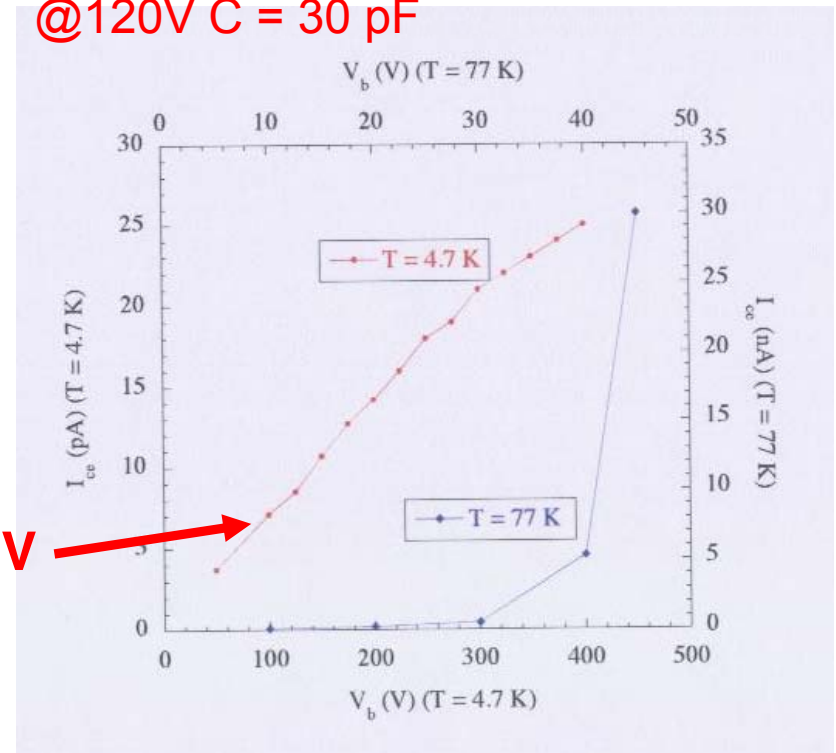

Cold Electronic Chain

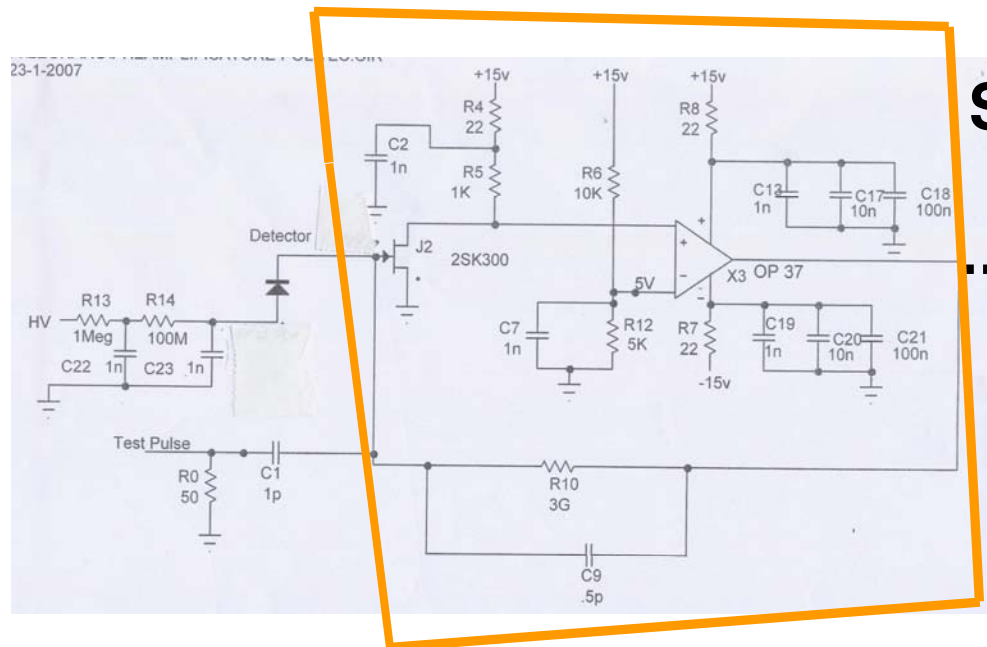

Si Detector at 4,7 Kelvin


To Silena Shaping Amplifier

EGG MCA

Charge Amplifier at 120 K


Electronic Noise = 1,9 KeV

FWHM (Right side) 60 Kev =4,2 KeV

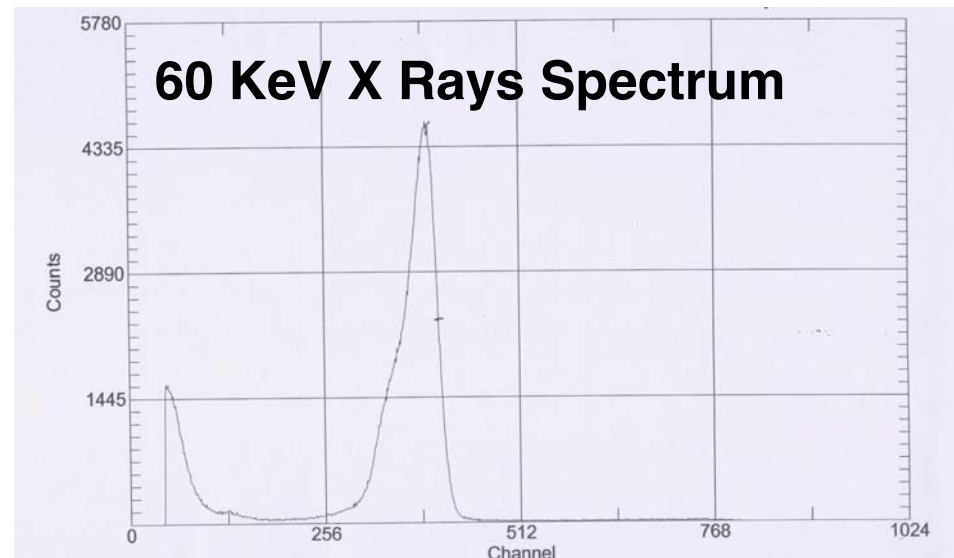

Peak Positions vs Vbias



----- 60 KeV + Cosmic Rays -----

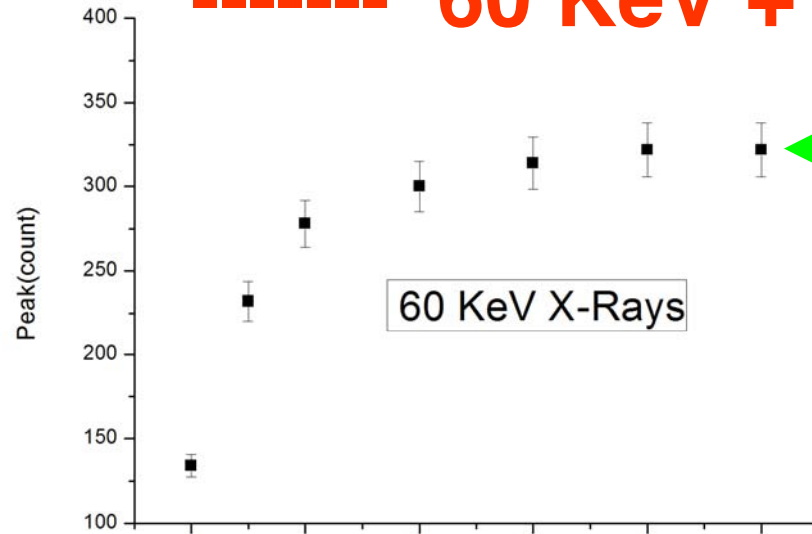

100 % Charge Collection Efficiency

PIN diode P-type 40 KOhm*cm

Plateau value at 90 - 100 V/cm

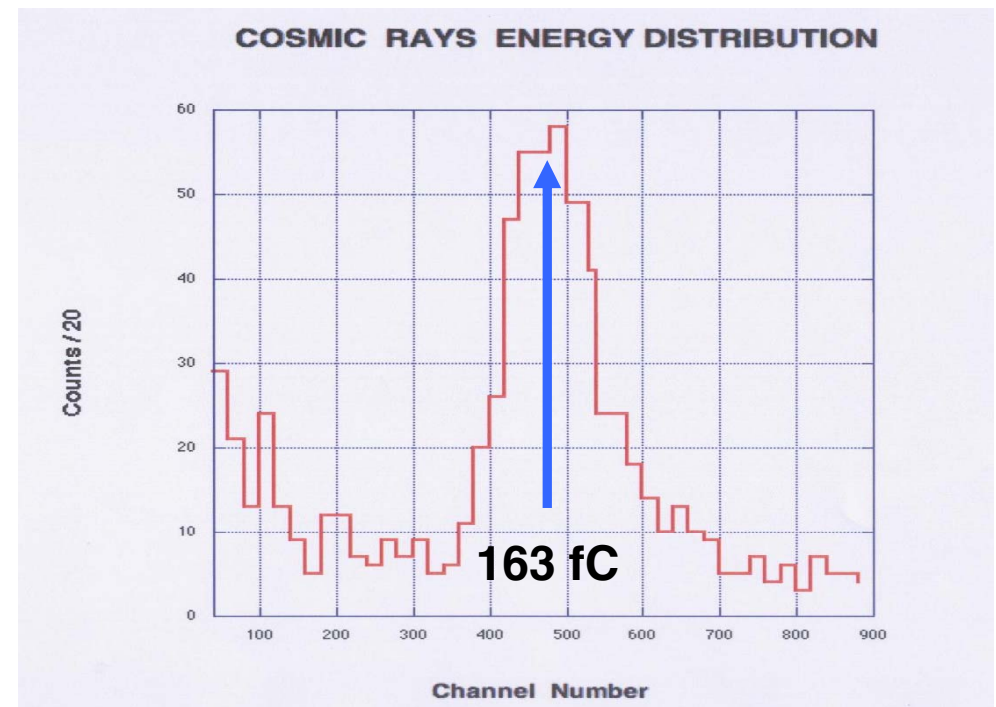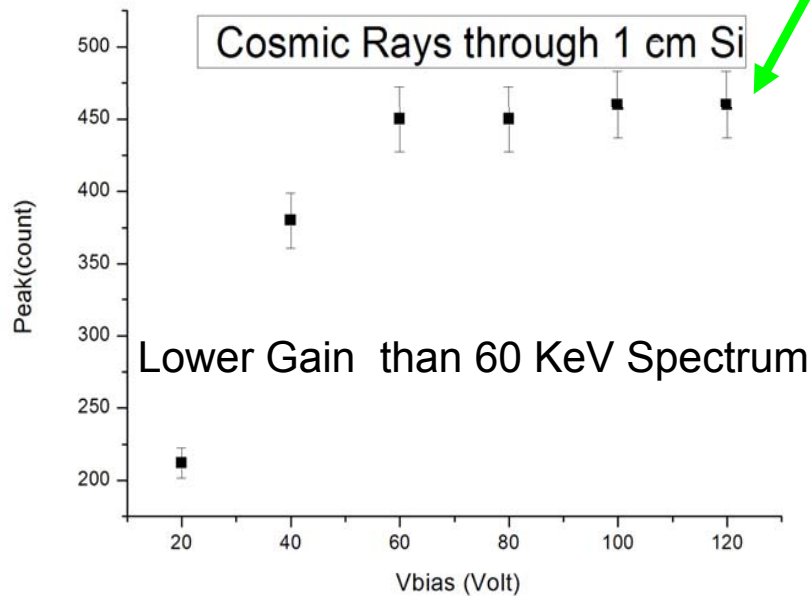

Collection efficiency vs Purity


E (V/cm) @ 100% charge coll.eff. VS **Resistivity (KOhm*cm @ RT)**

a) 1000	6,8	n - type	Wacker
b) 125	12	n – type	Wacker
c) 90	40	p – type	Top-Sil

- a) Braggio et al NIMA 568 (2006) 412
- b) Braggio et al NIMA accepted for pub.
- c) Braggio et this work

Charge Collection Efficiency Depends Critically from Purity Sample

Home Made Germanium Schotky Barrier Detector

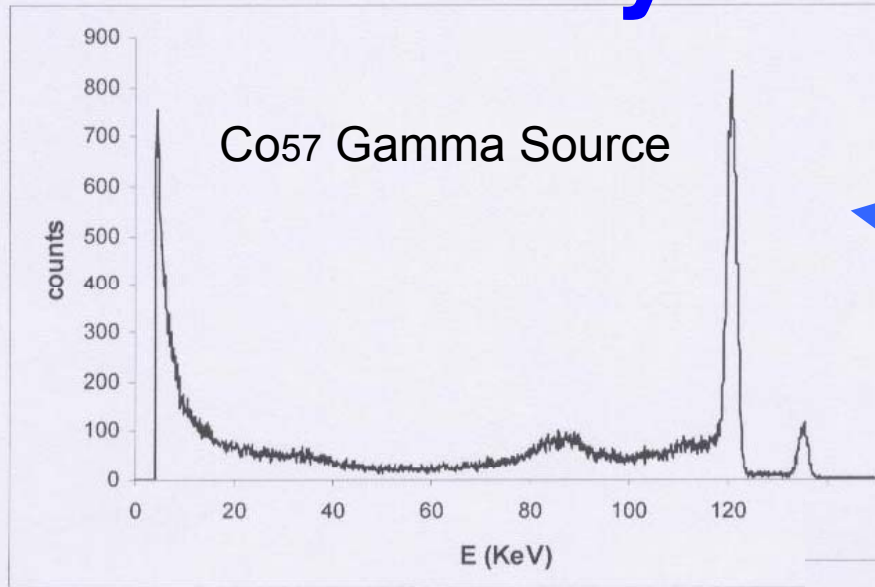

Cleaning procedure and electrode evaporation made in our lab.


1 Cm thick 3 cm Diam.

Vbias = 25 Volt

Germanium purity 2 OfM Better than Silicon

Germanium detector cost 50 times
Lower than same planar device

Energy Resolution comparable
With a standard Germanium

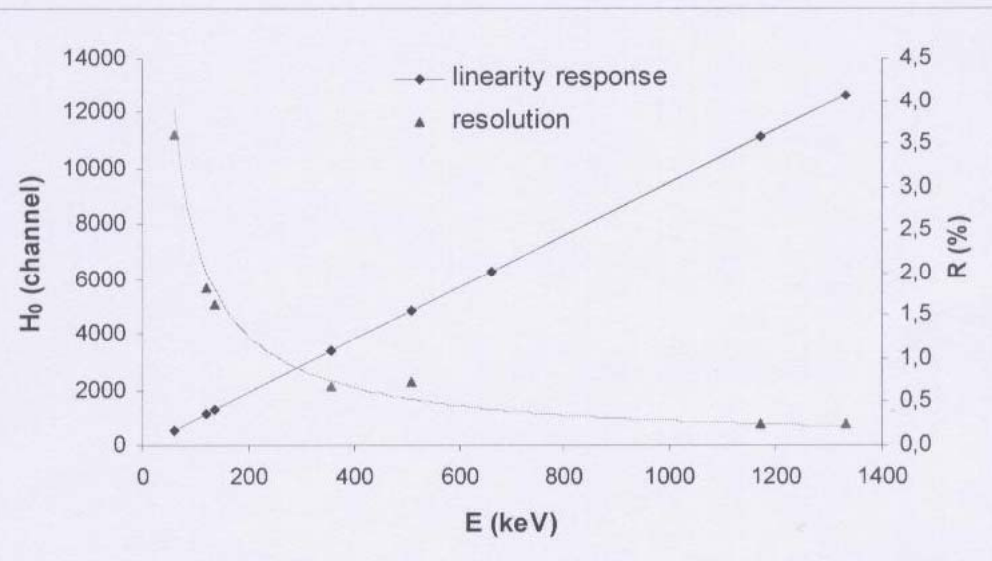

Conclusions

Perspectives at 4K

- 1) Freeze-out effect observed clearly through Capacity at LHe Temperature
- 2) Leakage Current less than 1 pA/cm^2
- 3) Low E Field 100% Charge Collection Efficiency depends on Samples Purity
- 4) Detectors Processing easy to arrange also with X-Y Strips
- 5) High Purity Samples can be processed to have 5 cm Thick Detectors TPC
- 6) Large electron mean free path decreases the Avalanche Electric Field
- 7) 4 Kelvin Cryocooler available nowadays No Liquid Helium handling